

TH Premium SR LI 19

Technical Data

Made for
iPhone | iPad | iPod



- 60 dB / 119 dB SPL
(2 ccm coupler)
- 70 dB / 129 dB SPL
(ear simulator)

TH Premium SR LI 19 | Technical Data

Type

Receiver



Output sound pressure level

OSPL 90 at 1.6 kHz

OSPL 90 (Peak)

HFA-OSPL 90

Gain

FOG at 1.6 kHz

FOG (Peak)

HFA-FOG

Reference test gain

Frequency, noise and directivity

Frequency range

Equivalent input noise

Total harmonic distortion at
500 / 800 / 1600 / 3200 Hz

Tinnitus therapy broadband

AI-DI

Inductive coil sensitivity

MASL (1 mA/m) at 1.6 kHz

HFA MASL (1 mA/m)

HFA SPLITS (left/right)

RSETS (left/right)

HFA SPLIV

Battery

Battery voltage

Battery current drain

Battery runtime (without streaming)

Battery runtime (incl. 5 h streaming)

IRIL IEC 60118-13:2016 Ed. 4.0

700-960 MHz (rating)

1400-2000 MHz (rating)

2000-2700 MHz (rating)

ANSI C63.19-2011

800-950 MHz (rating)

1600-2500 MHz (rating)

2 ccm coupler

Ear simulator

–

123 dB SPL

119 dB SPL

129 dB SPL

113 dB SPL

–

–

55 dB

60 dB

70 dB

50 dB

–

36 dB

48 dB

100 - 9000 Hz

100 - 10000 Hz

19 dB SPL

24 dB SPL

2 / 2 / 2 / 1 %

3 / 5 / 4 / - %

70 dB SPL

–

4.0 dB

–

–

–

–

–

–

–

–

–

–

1.25 V

1.4 mA

1.4 mA

up to 19 h

up to 16 h

user

user


user

M4

M4

Please find additional information to the values on page "Further Information"


TH Premium SR LI 19 | Fitting Range

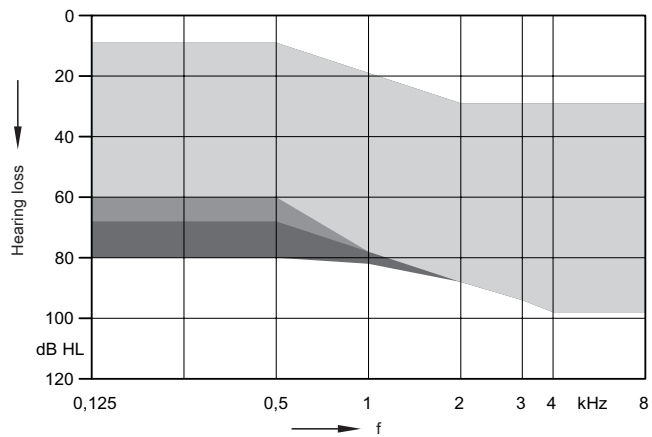

- Open Click Domes
- Double Click Domes
- Click Mold (no vent)

Receiver (Closed Click Dome) | Basic Data


2 ccm coupler

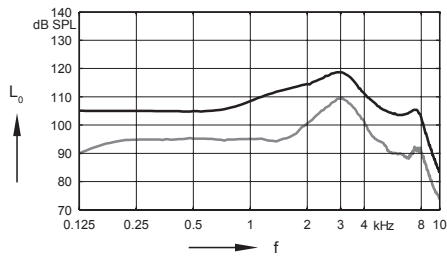

Max. Output sound pressure level
($L_1 = 90$ dB)


Full on gain
($L_1 = 50$ dB)

Ear simulator

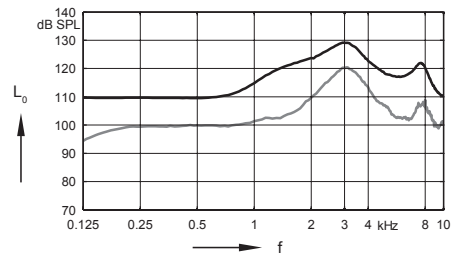


Max. Output sound pressure level
($L_1 = 90$ dB)

Full on gain
($L_1 = 50$ dB)

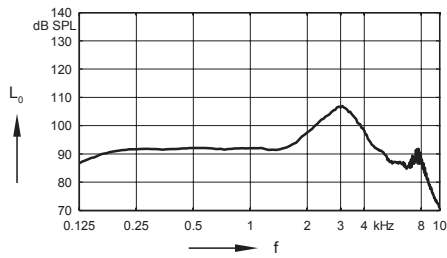

Frequency response
($L_1 = 60$ dB)

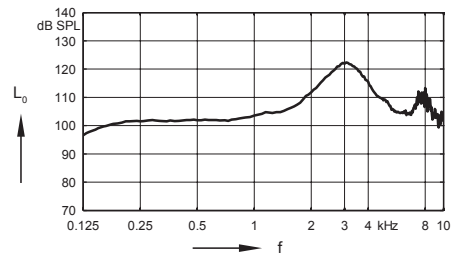

Basic acoustic response
($L_1 = 60$ dB)

TH Premium SR LI 19 | Features and Accessories

Audiology	
Own Voice Processing (OVP) ¹⁾	■■■■■
3D Classifier	■■■■■
Signal processing (channels) / Gain/MPO (handles)	48 / 20
Hearing programs	6
Sound Clarity	
HD Spatial	●
Extended dynamic range	●
Extended bandwidth	●
EchoShield	●
HD Music (presets)	3
eWindScreen binaural ¹⁾	●
eWindScreen®	●
Noise Management	
Speech and noise management (steps)	7
SoundSmoothing (steps)	3
Directional speech enhancement (steps)	3
Feedback cancellation	●
Speech Quality	
Directionality	
Narrow Directionality ¹⁾	●
Spatial SpeechFocus ¹⁾	●
SpeechFocus	●
TwinPhone ¹⁾	●
Frequency compression	●
Direct Streaming	
Made for iPhone	●
Adaptive Streaming Volume ²⁾	●
Tinnitus	
Notched Noise Therapy	●
Tinnitus therapy	●
Fitting	
Smart Optimizer and Data Logging	●
Acclimatization manager	●
Performance Guide	●
Insitugram	●
Learning (classes)	6
TeleCare	
TeleCare 3.0	●

¹⁾ req. bilateral fitting

²⁾ streaming only

● available ■■■■■ highest feature performance — not available

TH Premium SR LI 19 | Features and Accessories

Style specific features	
Ingress Protection Rating	IP68
Charging contacts	●
Battery Size	—
Battery door on/off function	—
Nanocoated housing	—
e2e wireless® 3.0	●
User controls coupling via e2e	●
Wireless programming	●
Instrument configurations	
Flat cover	—
Rotary volume control	—
Push button	—
Rocker switch	—
Color conversion kit	—
Battery door – integrated telecoil	—
Battery door – child lock	—
Small earhook	—
Programming accessories	
ConnexxAir / ConnexxLink®	— / —
NoahLink wireless	●
Programming adapter / cable	—
Accessories	
miniPocket	○
TH CROS RIC 19	—
TH CROS LI 19	—
StreamLine TV	○
StreamLine Mic	○
Inductive Charger	mandatory
Apps	
myControl App	○
touchControl App	○

● available ○ optional — not available

Further Information

Abbreviations

The following abbreviations are used in this datasheet:

OSPL	Output Sound Pressure Level
HFA	High Frequency Average
FOG	Full On Gain
MASL	Magneto Acoustical Sensitivity Level
SPLITS	Coupler SPL for an Inductive Telephone Simulator
RSETS	Relative Equivalent Telephone Sensitivity
SPLIV	SPL In a Vertical magnetic field
AI-DI	Articulation Index - Directivity Index
IRIL	Input Related Interference Level
RTF	Reference Test Frequency

Standards and additional information

- ▶ All measurements with the 2 ccm coupler were performed according to ANSI S3.22-2014 and IEC 60118-0:2015 if applicable.
- ▶ All measurements with an ear simulator were performed according to IEC 118-0/A1:1994 and to DIN 45605 (frequency range) if applicable.
- ▶ Curves and figures representing FOG are measured with 20 dB reduction and 70 dB SPL input level.
- ▶ Figures representing Equivalent Input Noise incorporate a moderate expansion.
- ▶ Tinnitus therapy measurement conditions: all tinnitus single frequency sliders in max position, master volume slider in default position (0 dB) and local volume control in default position.
- ▶ Inductive coil sensitivity values, inductive response curves and T ratings apply for instruments with telecoil only.
- ▶ The following acoustic connections / ear pieces were used:
 - Receiver Unit: Closed Click Dome
- ▶ The current consumption is measured in reference test setting (RTS) according to the applicable standards. Due to the settling behaviour of hearing instruments supporting RF (radio frequency), the battery current is measured 3 minutes after turning on (note: no pairing).
- ▶ The battery runtime is based on first fit settings using 60% of the fitting range and an ISTS (International Speech Test Signal) input signal at 65 dB SPL (note: pairing established). The actual battery runtime is determined by battery quality, hearing loss, sound environment, usage and activated feature set. Regarding RF usage (Bluetooth streaming) two different conditions are considered.

Special note for instruments with built-in lithium-ion rechargeable battery

- ▶ The runtime of all lithium-ion rechargeable batteries reduces over time. The estimates are based on fresh lithium-ion rechargeable battery capacity. Under normal operating conditions, the battery will retain up to 80% of its initial capacity after 2 years of use. Please note that battery performance will vary depending on individual usage patterns and environmental conditions.



“Made for iPod”, “Made for iPhone”, and “Made for iPad” mean that an electronic accessory has been designed to connect specifically to iPod, iPhone, or iPad, respectively, and has been certified by the developer to meet Apple performance standards. Apple is not responsible for the operation of this device or its compliance with safety and regulatory standards. Please note that the use of this accessory with iPod, iPhone, or iPad may affect wireless performance.

The information in this document contains general descriptions of the technical options available, which do not always have to be present in individual cases and are subject to change without prior notice. The required features should therefore be specified in each individual case at the time of conclusion of the respective contract.

Manufactured for

TruHearing Inc.
12936 S. Frontrunner Blvd
Draper, UT 84020
United States

Order No. 03999-99T1-7600 SI/19327A-20
© 08.2019, TruHearing Inc.
All rights reserved

www.truhearing.com



Warning

Choking hazard posed by small parts.

- ▶ This instrument is not intended for the fitting of infants, children under 3 years and persons of mental incapacity.



Warning

Instrument has an output sound pressure level of 132 dB SPL or more.

Risk of impairing the residual hearing of the user.

- ▶ Take special care when fitting this instrument.